

ISOSAN

100% STEAM SANITISATION



PROMED
MEDICAL COMPANY



ISOSAN

THE PROFESSIONAL SOLUTION AGAINST VIRUSES, BACTERIA AND ALLERGENS

ISOSAN by Promed is designed to sanitise equipment, surfaces and whole rooms using a combination of dry saturated steam and disinfectant*. Able to reach a steam temperature of 165°C, **it kills 99.999% of viruses, germs, bacteria, fungi and spores**. The equipment surrounds surfaces with a fast-drying fog. Made entirely in Italy, with a steel construction and copper boiler with external heating element, it saves energy during use and lasts longer as the heating element is not immersed in water. Reaching a power output of 3kW, it reaches the work temperature more quickly (in just 7 minutes) and delivers 6 bar pressure. It uses tap water, which is added to the non-pressurised tank enabling constant sanitisation cycles with no interruptions.



PM-BF009PL3

STEAM COMBINED WITH DISINFECTANT*

By fitting the "steam fogger" accessory, steam combines with detergent to produce sanitising fog. In contact with surfaces, the fog cools, condenses into disinfectant micro-droplets, and creates a protective sanitising film that guarantees the **immediate, deep-down and long-lasting elimination of viruses, bacteria and other pathogens**.


ACCESSORIES PROVIDED



Steam handle 5 m



D28 nylon bristle brush for steam nozzle

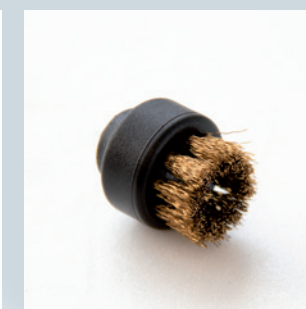

D28 brass bristle brush for curved steam nozzle




Steam nozzle



Curved steam nozzle



Steam fogger

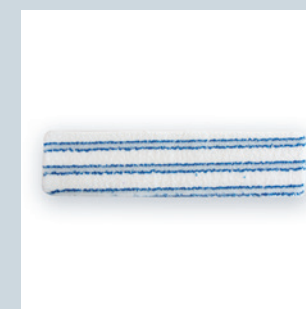

PM-RIP0811
Blue striped microfibre cloth 400x100 hard surfaces for RIP0801



PM-RIP0801
Professional steam brush



PM-ISO102310

ISOBACTY

SANITISER 70% ALCOHOL + HYDROGEN PEROXIDE


Ultra-rapid, ready to use and high efficacy alcoholic solution. The high concentration of pharmaceutical grade ethanol, in tandem with hydrogen peroxide (active oxygen), guarantees deep-down, rapid hygiene on all surfaces without leaving residue.

ISOSAN

100% STEAM SANITISATION

TECHNICAL DATA

POWER SUPPLY:	220V/240V 50/60Hz
BOILER CAPACITY:	2,8 L
BOILER POWER:	3000W
STEAM PRESSURE:	6 Bar
STEAM TEMP. IN BOILER:	165°C
STEAM OUTPUT:	60 g/min
WATER REFILL TANK:	3 L
DIMENSIONS:	28,5x38xh48 cm
WEIGHT WITHOUT ACCESSORIES:	11 kg
WEIGHT WITH ACCESSORIES:	13 kg
ADJUSTABLE STEAM SUPPLY	✓
BOILER REFILL	AUTOMATIC
WATER OUT AUDIBLE ALARM	✓

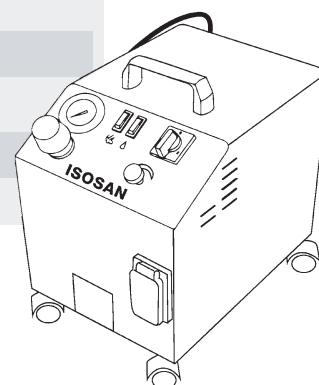

EFFICIENCY, SAFETY AND NOTES

Reliability and **safety** guaranteed by a **copper boiler**. Copper is flame-proof and anti-limescale. **Heating elements** are external, so **never** come in contact with water.

1 kW power generates a steam flow rate of 1.2 kg/h (20 g/min). The effectiveness of steam is proportional to flow rate, not pressure. High pressure teamed with low flow rate is ineffective.

MADE IN ITALY

PROMED

M E D I C A L C O M P A N Y

Promed srl - Via Rossini, 26
40055 Villanova di Castenaso (BO), Italy
Tel. +39 051 0390446 - sales@promeditaly.it
www.promeditaly.it